

Container Filter Konteyner Filtresi

@x
AXIS

Product Code/Ürün Kodu
AXCP205



100 pcs/adet

compatible with ISO 11140-1:2017
ISO 11140-1:2017 ile uyumludur.

Manufacturer / Üretici
KORDON TIP SAĞLIK ARAC GEREÇLERİ
MUH.PROJE.ITH.SAN.VE TIC.LTD.STİ
10006/1 Sok. No:43 35620 A.O.S.B CİGLİ İZMİR TURKEY
+90 232 348 60 60 / +90 232 348 60 90
www.kordontip.com - info@kordontip.com



Container Filter

for steam sterilization processes

Product Description

Container filters for Steam sterilization processes.

Initial color: pink

Final color: siyah

Container labels are composed of a identification card with an chemical indicator. The area provides the space for identification of sterilizer no, sterilization date, due date an name. When label exposed to steam the pink indicator turns black.

Application

The product is used to identify sterilized materials in the sterilizer. It is placed on the container according to the regulations of the institution or hospital. The yellow indicator on the label turns grey when exposed to steam. It is placed on the container in order to ensure the exposure and sterilization safety of the container after sterilization.

Storage

Indicators should be stored in dry places, protect from light at temperatures between 15 - 30 °C, 35 - 60 % relative humidity. Do not storage close to sterilizing agents. Do not use indicators after their expiration date. Do not reuse.

Not processed



Unsatisfactory



Processed
and successful



Konteyner Filtresi

buhar sterilizasyon işlemlerinin izlenmesi için

Ürün Açıklaması

Buharlı sterilizasyon işlemleri için konteyner filtresi.

İlk renk: pembe

Son renk: siyah

Konteyner etiketleri, kimyasal indikatörlü bir kimlik kartından oluşur. Alan, sterilizatör numarasının, sterilizasyon tarihinin, vade tarihinin ve adın tanımlanması için alan sağlar. Etiket buhara maruz kaldığında pembe indikatör siyaha döner.

Uygulama

Ürün sterilizatörde steril edilmiş malzemelerin tanımlanması için kullanılmaktadır. Kurumun veya hastanenin düzenlemelerine göre konteyner üstüne yerleştirilir. Etiket üzerindeki sarı indikatör buhara maruz kaldığında griye döner. Konteynerin sterilizasyondan sonra maruziyetini ve sterilizasyon güvenliğini sağlamak amacıyla konteyner üzerine yerleştirilir.


Depolama

İndikatörler kuru, ışıktan korunan, 15 - 30 °C sıcaklıkta, %35 - 60 bağıl nemde saklanmalıdır. ıslatmayın. Sterilize edici maddelere yakın saklamayın. Son kullanma tarihi geçen göstergeleri kullanmayın. Tekrar kullanmayın.

İşlem öncesi



Hatalı

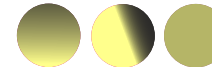

İşlem görmüş
ve başarılı



ENG

TR